


Fatty Acids (Free) Rtx[®]-200

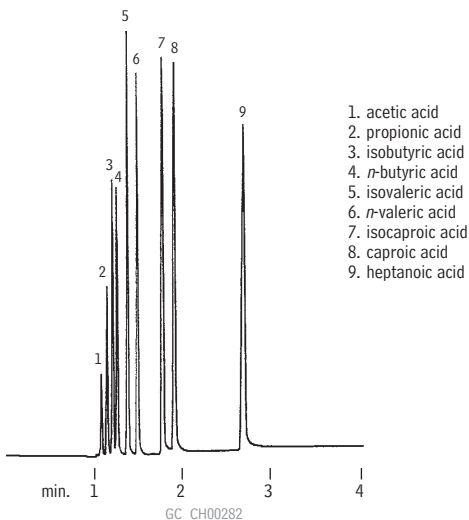

30m, 0.25mm ID, 0.25 μ m Rtx[®]-200 (cat.# 15023)
Sample: 0.8 μ L split injection of a free fatty acid standard.
Conc.: approximately 10 to 20ng/ μ L.
Oven temp.: 90°C
Inj. & det. temp.: 250°C
Carrier gas: hydrogen
Linear velocity: 40cm/sec. (flow rate: 1.4cc/min.)
FID sensitivity: 4 x 10⁻¹¹ AFS
Split vent: 40cc/min.