


Organic Acids Free Fatty Acids Stabilwax®-DA

(direct injection)

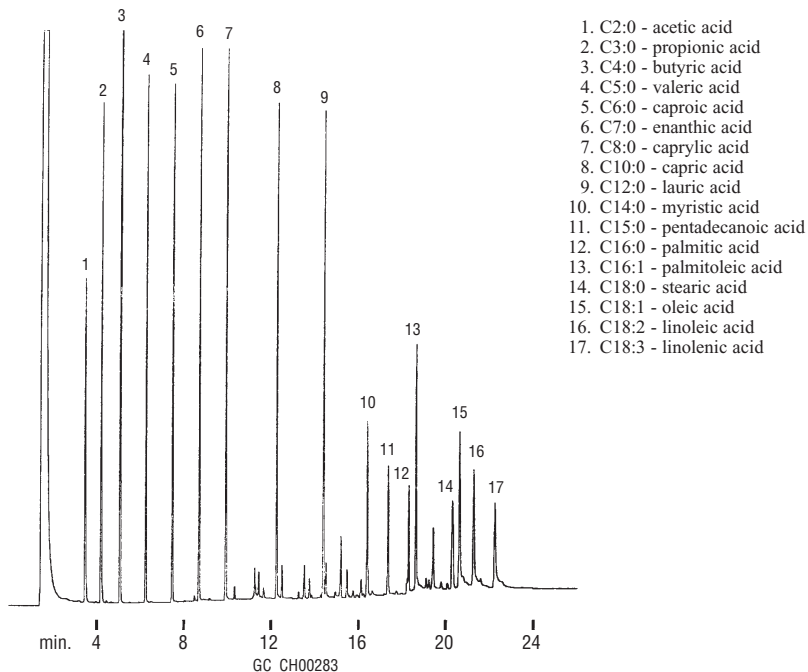

30m, 0.53mm ID, 0.25µm Stabilwax®-DA (cat.# 11025)
0.5µL direct injection of a 5mg/mL standard.

Oven temp.: 100°C (hold 2 min.) to 250°C @ 8°C/min.
Inj. & det. temp.: 280°C
Carrier gas: helium
Linear velocity: 40cm/sec. (flow rate: 5.2cc/min.)
FID sensitivity: 8×10^{-11} AFS