


Chlorofluorocarbons

Rtx®-1

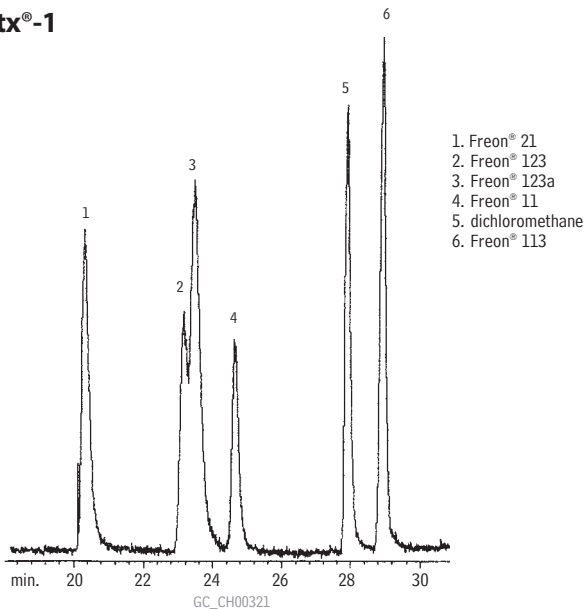

105m, 0.53mm ID, 3.0 μ m Rtx®-1 (cat.# 10189)
1.0 μ L splitless injection of a 500 μ g/mL mix of Freon® compounds

Splitless hold: 0.5 min.
Oven temp.: 35°C
Inj. temp.: 220°C
Carrier gas: helium
Linear velocity: 20cm/sec. set @ 35°C
MSD: TIC mode
Scan rate: 1.5 scans/sec.

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com