


## Solvents Rtx®-Volatiles

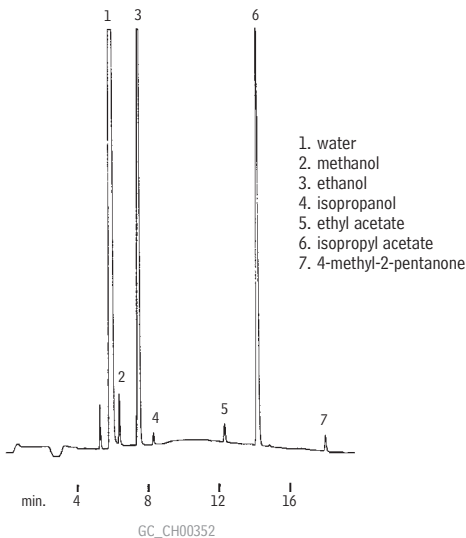

60m, 0.53mm ID, 2.0 $\mu$ m Rtx®-Volatiles (cat.# 10905)  
1.0 $\mu$ L split injection of neat solvents

Oven temp.: 40°C (hold 2 min.) to 200°C @  
5°C/min. (hold 10 min.)  
Inj. & det. temp.: 250°C  
Carrier gas: helium  
Linear velocity: 20cm/sec. (flow rate: 2.6cc/min.)  
TCD sensitivity: 4 x 10<sup>-11</sup> AFS  
TCD current: 100mA  
Split ratio: 40:1

Restek Corporation 110 Benner Circle Bellefonte, PA 16823  
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com