


# Cleaning Solvents (All-Purpose Cleaner) Rtx®-VMS

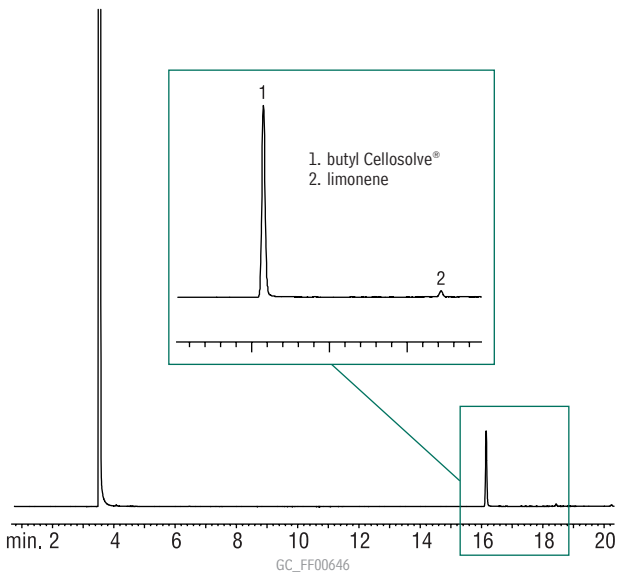

Rtx®-VMS, 60m, 0.25mm ID, 1.4 $\mu$ m (cat.# 19916)  
Sample: approx. 5% cleaner in methanol  
Inj.: 1.0 $\mu$ L split (split ratio 100:1), 4mm inlet liner (cat.# 20814)  
Inj. temp.: 250°C  
Carrier gas: hydrogen, constant flow  
Flow rate: 1.5 mL/min.  
Oven temp.: 60°C (1 min. hold) to 138°C @ 4°C/min. to 145°C @ 2°C/min. to 210°C  
@ 15°C/min. to 230°C @ 5°C/min.  
Det.: FID @ 250°C

Restek Corporation 110 Benner Circle Bellefonte, PA 16823  
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • www.restek.com