


# Underivatized Ibuprofen Rt- $\beta$ DEXsm™

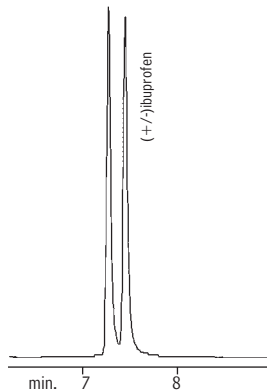

GC\_PH00242

30m, 0.32mm ID, 0.25 $\mu$ m Rt- $\beta$ DEXsm™ (cat.# 13104)

On-column injection

Conc.: 125ng each enantiomer

Oven temp.: 175°C to 200°C @ 2°C/min.

Inj./det. temp.: 200°C/230°C

Carrier gas: helium

Linear velocity: 60cm/sec.

Det.: GC-FID

Split ratio: 13:1 using cup splitter liner (cat.# 20709)

Restek Corporation 110 Benner Circle Bellefonte, PA 16823  
814-353-1300 • 800-356-1688 • Fax: 814-353-1309 • [www.restek.com](http://www.restek.com)