


Underivatized Steroids, Anabolic Rtx®-5

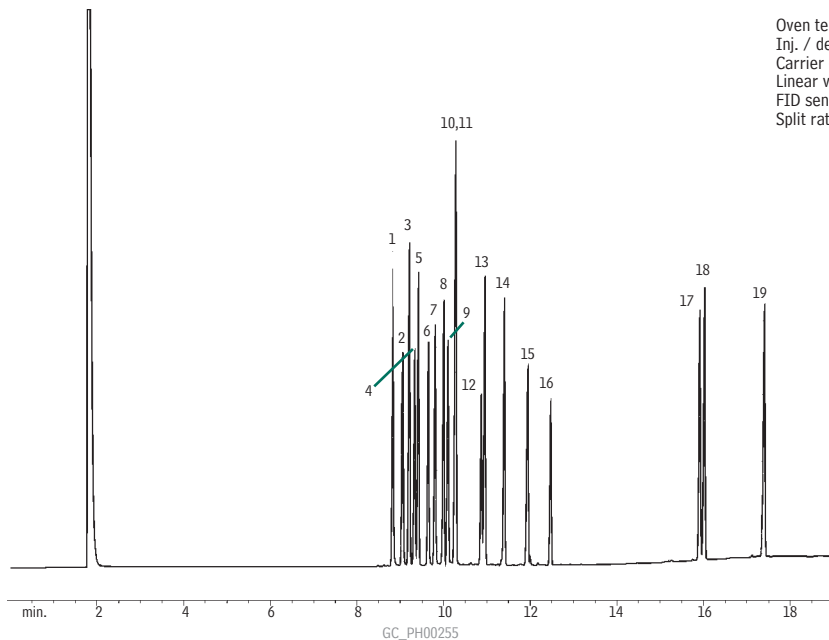

30m, 0.25mm ID, 0.10 μ m Rtx®-5 (cat.# 10208)
0.5 μ L split injection of anabolic steroids
Conc.: 1000ng/ μ L

Oven temp.: 180°C to 340°C @ 10°C/min. (hold 3 min.)
Inj. / det. temp.: 280°C / 340°C
Carrier gas: helium
Linear velocity: 35cm/sec. set @ 180°C
FID sensitivity: 1.28 x 10¹⁰ AFS
Split ratio: 50:1

1. 5-androstene-3 β ,17 β -diol
2. 17 α -methyl-5-androstene-3 β ,17 β -diol
3. 5 α -androstan-17 β -ol-3-one
4. 19-nortestosterone
5. 17 α -methylandrostan-17 β -ol-3-one
6. mesterolone
7. testosterone
8. 17 α -methyltestosterone
9. 1-dehydrotestosterone
10. 1-dehydro-17 α -methyltestosterone
11. bolasterone
12. oxymethalone
13. 19-nortestosterone-17-propionate
14. testosterone propionate
15. fluoxymesterone
16. 4-chlorotestosterone-17-acetate
17. testosterone-17 β -cypionate
18. 1-dehydrotestosterone benzoate
19. 1-dehydrotestosterone undecylenate