


Underivatized Phenothiazines (Basic Drugs) Rtx[®]-5

1. phenothiazine
2. promethazine
3. trimeprazine
4. promazine
5. chlorprothixene
6. chlorpromazine
7. trifluoperazine
8. acetopromazine
9. prochlorperazine
10. thioridazine

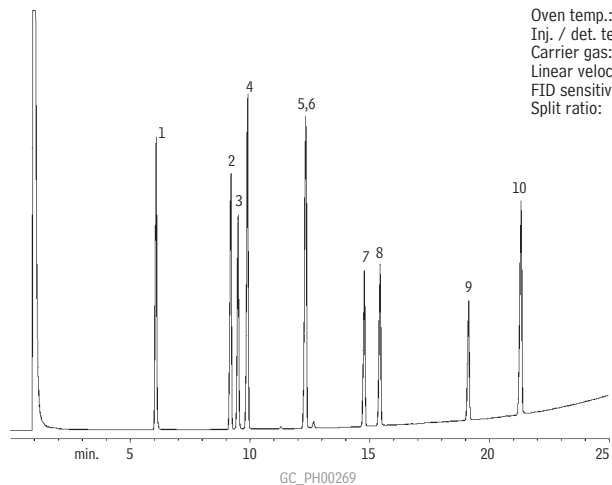

15m, 0.32mm ID, 0.50 μ m Rtx[®]-5 (cat.# 10236)
1.0 μ L split injection of phenothiazines
Conc.: 2000ng/ μ L

Oven temp.: 200°C to 325°C @ 5°C/min.
Inj. / det. temp.: 250°C / 315°C
Carrier gas: helium
Linear velocity: 25cm/sec. set @ 200°C
FID sensitivity: 2.56 x 10⁻¹⁰ AFS
Split ratio: 30:1